

Windows NT4.0 Printer Driver on Windows2000 (Parallel Connection)

*This explanation is the case of the UP-D890 printer driver installation.

(Installation)


1. Connect the printer to your PC.
2. Power on the printer.
3. Start up the Windows2000.
4. The following dialog is displayed during startup. (If not, please install from [Add printer] in [Printers] folder.)




5. Click on [Next].




6. Choose [Display a list of the known drivers for this device so that I can choose a specific driver] and click on [Next].

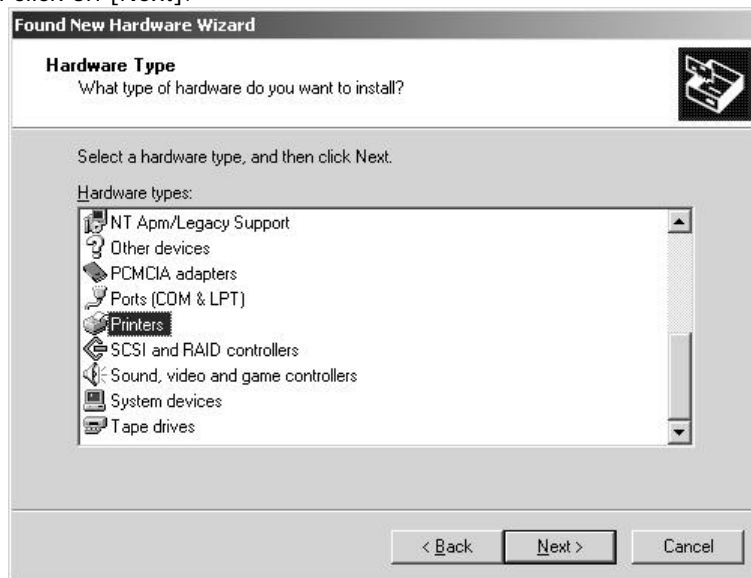

7. Choose [Printers] and click on [Next].

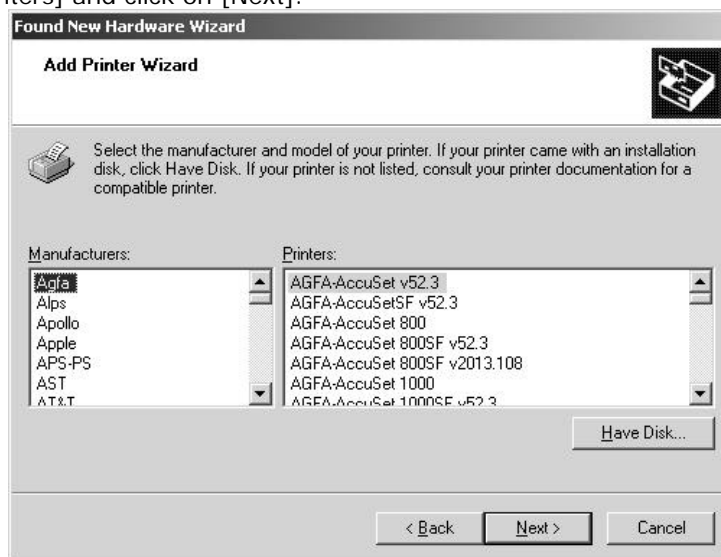


8. Click on [Have Disk...].

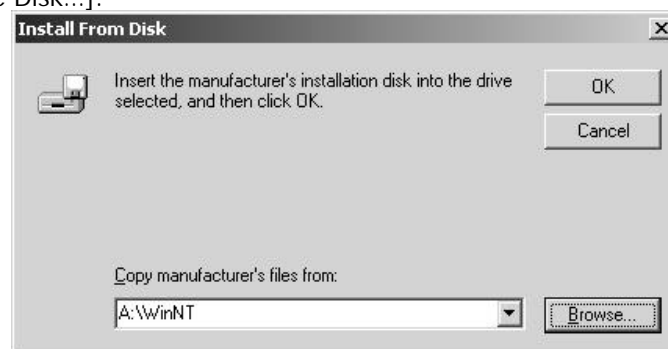

9. Specify the folder in which the printer driver is and click on [OK].




10. Choose [Sony UP-D890] and click on [Next].

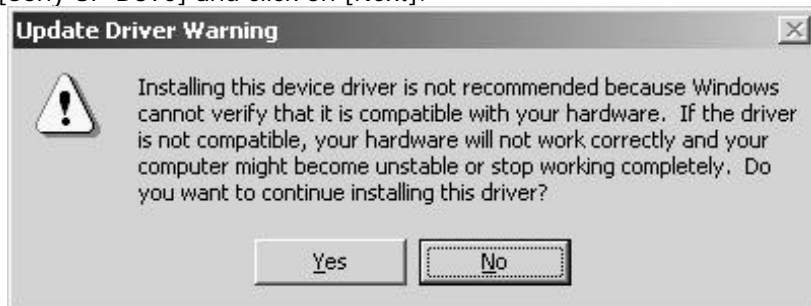

11. Click on [Yes].



12. Click on [Next].

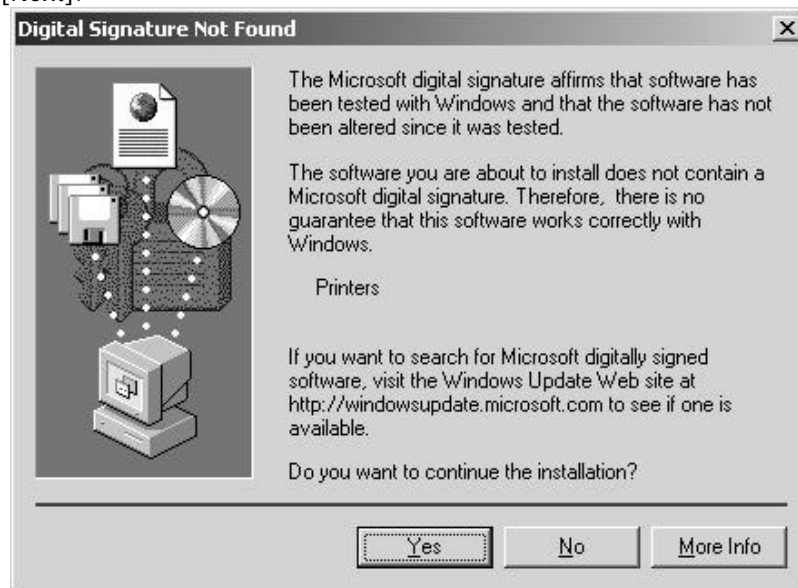


13. Click on [Yes].



14. Click on [Finish]. This completes the installation.

(Port setting)

Bi-directional communication does not work correctly. So, please do not click [Enable bidirectional support] in [Ports] in [Properties] of the printer.

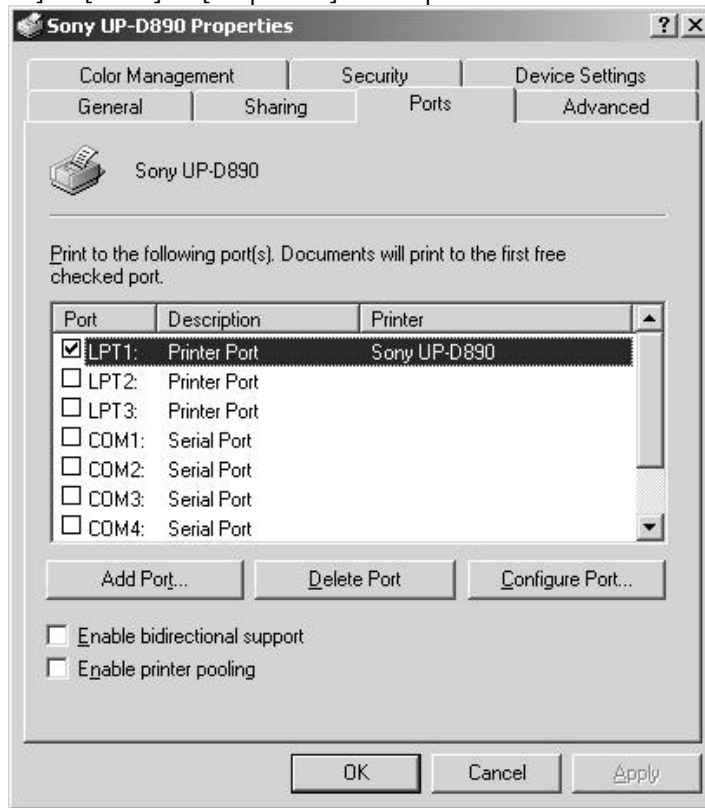

In case you OFF the bidirectional communication as above, there are following limitations.

- Detailed status error is not displayed.
- Even if you cancel the jobs in spooler, printer is not reset. To start printing again, please restart the printer.

(About ECP mode)

Even if you set the BIOS of your PC to ECP, it does not work in ECP mode. Now the cause is under investigation.